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— ENGLISH —

IMPORTANT SAFETY WARNINGS

- Do not allow any foreign objects to enter the clearance between the massage head and the casing.
- Do not allow the product to be exposed to water to avoid any possible accidents that may result from a short circuit.
- If the mains cable or plug is damaged, do not disassemble and attempt to repair it by yourself. Please contact your distributor.
- Ensure that the mains power is suitable for this product.
- Please disconnect the product from the mains supply after use or before cleaning to avoid damage to the equipment or personal injury.
- Please use this product according to this Manual.
- Do not use any spare parts and attachments other than those recommended.
- Do not use this product if its casing or mains cable is damaged.
- Do not use this product if you are under the influence of alcohol or feeling unwell.
- Do not use this product with too much force to avoid injury.
- Do not keep this product on for more than 15 minutes to avoid overheating.
- Do not use this product within one hour of eating.
- Please use this product carefully according to your doctor's advice if you have a neck or back injury or if you have had an operation recently.
- This appliance contains batteries that are non-replaceable.
- WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
- Adapter: Xiamen Invcy Electronics Co., Ltd. IVP0900-2000G

WORKING ENVIRONMENT

- Do not use this product in a highly humid environment such as a bathroom.
- Do not use this product immediately after there has been a sharp change in the environmental temperature.
- Do not use this product in an environment where there is a lot of dust or any corrosive gases.
- Do not cover this item by a blanket or pillow, etc. during use.

PEOPLE THAT SHOULD NOT USE THIS PRODUCT

- People with contagious skin conditions.
- Patients with osteoporosis.
- Patients with heart diseases and with electronic medical instruments like a heart pacemaker.
- People with fever.

- Women that are pregnant.
- Those who are injured or suffering from any body surface disease.
The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Children under 12-years old.
- This device is not suitable for use by persons (including children) with reduced physical, sensory or mental capabilities or who have limited experience and knowledge, unless they are under the supervision of, or have been trained to use the device, by a person responsible for their safety.

SAFETY AND MAINTENANCE ISSUES

- Precautions must be taken to prevent children playing with the product.
- Those who to need rest according to a doctor's advice or those who do not feel well.
- Do not use this product when you are wet.
- Do not use this product on a pet.

SAFETY PRECAUTIONS

- Do not plug or unplug this product when your hands are wet.
- Do not pull the cable when unplugging. Do not plug or unplug carelessly.
- This is not a toy. Keep away from children.
- Do not damage the cable or change the circuit of the product.
- Please stop using this product in case of a power cut in order to avoid any body injury upon sudden power resumption.
- Please stop using this product and consult your local supplier if you notice anything abnormal with the product during use.
- Please stop using this product and consult your doctor for advice
- If the supply cord is damaged, the manufacturer, its service agent or similar qualified persons must replace it, in order to avoid a hazard.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- When pulling out power plug, Please Hold the plug pulled back, rather than directly draw the power cord.
- When the whole machine does not working, please immediately stop using it to pull out power plug.
- Do not use this item in the following locations: Bathroom and kitchen, high temperature lights, damp, high voltage and the source of fire and water.
- Keep away from sharp objects, the machine cannot be used when its textile cover is damaged.
- Do not force it with high pressure or put body pressure on the massage heads. Do not touch with hard objects.

PRODUCT FEATURES

Armrest

Heating

- This massager is designed to fit your neck and shoulder, but is also useful to massage your back, thigh and legs.
- Its shape and arm in arm design give you the choice to how strong or gentle massage you want.
- You can select clock or anti-clockwise circles of the massage heads.
- The warm function promote's the blood circulation.


1. Neck & Shoulder

2. Abdominal Massage

3. Waist Massage

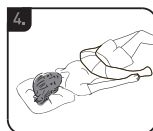4. Thigh Massage

CONTROL PANEL

Heat button

Power button

Kneading direction

Charging light - Status
Red - Charging, Green - Fully charged

HEAT

Turn on or off the heat function.

ON/OFF


Start and stop the massage. After start, the massage automatically stop after 15 minutes.

DIRECTION

Change the direction of the massage heads.

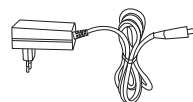Plug the power socket to your power outlet.
Plug the charging plug to the massagers charging outlet.

LI-ION BATTERY WARNINGS

This device contains a Lithium-ion battery. The battery must be recycled or disposed of properly. Contact your local waste management office for information on battery recycling or disposal.

IMPORTANT: HANDLING OF RECHARGEABLE BATTERIES

- If you come in contact with fluid from the battery cell, rinse the affected areas with water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, open or crush the rechargeable batteries.
- Only use the attached charger.
- Batteries must be fully charged before first time use.
- Correct charging must be observed at all times.
- Fully charge the battery prior to initial use

CHARGING

Fully charge the massager before first time use.

1. Plug the AC ADAPTER into the DC IN JACK. Plug the AC ADAPTER into an electrical outlet.
2. The CHARGING LED INDICATOR will turn red to indicate charging.
3. The CHARGING LED INDICATOR will turn green when the massager is fully charged.
4. Approximate charging time is 150 minutes.

CLEANING AND MAINTENANCE

1. After use, let the massager cool down before you start cleaning.
2. After use, clean with damp cloth to remove dust and dirt from the surface.
3. Do not use corrosive cleaning solvent.
4. Only use dry cloth to clean control panel.
5. Store the massager in dry conditions.

TROUBLESHOOTING:

Problem	Cause	Solution
The item do not work	Battery need to be charged	Plug in charging adaptor and charge for 2 hours.
	The product has just finished 15 minutes of massage, and automatically turn off. The product has been overloaded. Too much weight on the massage heads.	Press on/off button to start again
		Turn off. Remove the weight. Unplug the power cord, wait until cool down. Turn on again.
The item is too hot	The heating function is on	Turn off the heating button
If none of this tips solve the problem, turn off the item for 60 minutes and make sure the battery is fully charged.		

CAUTIONS

Please do not use this product in dusty, hot or wet places.

This product must be disposed separately from other household waste as they can be harmful to the environment. Contact your local dealer or authority for more information.

ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE)

This symbol on the product or in the instructions means that the electrical and electronic equipment according to the European directive (WEEE) must be disposed of separately from other household waste, as they can be harmful to the environment. The countries in the EU have waste collection systems. Contact your local dealer of authority for more information.

CONTACT INFORMATION:

F&H Group A/S
Gl. Skivevej 70
DK-8800 Viborg
DENMARK