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F&H Group A/S
GL Skivevej 70
DK-8800 Viborg
DENMARK


SAFETY INSTRUCTIONS

When using electrical products, basic safety precautions should always be followed, including the following:

- The appliance is only intended for domestic/private use, not for commercial use.
- The appliance is only to be operated with the power adapter unit supplied, and with the mains voltage indicated on the power adapter unit.
- Never use the appliance if damaged or if any of its accessories are damaged.
- Never use the appliance during a thunderstorm.
- In the event of defect or malfunction, switch off the appliance immediately and disconnect from the power adapter and the mains outlet.
- Never pull the power cable when disconnecting from the power adapter or the mains outlet. Always pull the plug.
- Never hold or carry the appliance by the power cable.
- Keep cables away from hot surfaces.
- Always unplug the appliance from the electrical outlet after using and before maintenance. Turn all controls to the "OFF" position, and remove the plug from the outlet. - Never leave the appliance unattended when plugged in. Unplug from the mains outlet when not in use and before changing or taking off parts or attachments.
- The appliance can be used by persons from 12 years of age and above and persons with impaired physical, sensory or mental capabilities or lack of experience and knowledge if supervised or instructed in safe use of the appliance and they understand the hazards involved.
- Children must not play with the appliance.
- Do not reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug immediately. Keep dry- do not operate in wet or moist conditions.
- Do not insert fingers, pins, needles or sharp objects.
- Do not open the appliance, even using the zip fastener (fitted for reasons relating to manufacture.)
- Make sure all moving parts can move freely at all times.
- Do not sit or stand on the moving parts with all your weight as it may damage the appliance.
- A damaged cable must be replaced by authorized service personnel only.
- Never use under covers, such as pillows or blankets.
- Never use close to any flammable materials.
- Keep the appliance away from high temperatures.
- This appliance contains batteries that are non-replaceable.
- WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
- Adapter: Xiamen Invcy Electronics Co.,Ltd. IVP0900-2000G

PARTS

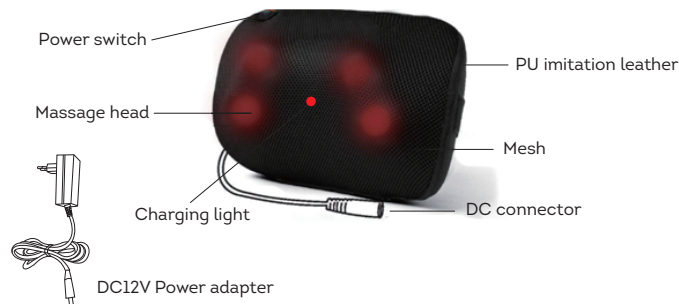LI-ION BATTERY WARNINGS

This device contains a Lithium-ion battery. The battery must be recycled or disposed of properly. Contact your local waste management office for information on battery recycling or disposal.

SPECIFICATIONS

Adapter: Input: 100~240V ~ 50/60Hz Output DC 9V 2A
Battery Type: Li-ion 2200mAh 7.4V
Battery Life: Approx. 120 min. (without heat) and approx. 70 min. (with heat)
Battery Charge Time: Approx. 150 min.

CHARGING

Fully charge the massager before first use.

1. Plug the AC ADAPTER into the DC IN JACK. Plug the AC ADAPTER into an electrical outlet.
2. The CHARGING LED INDICATOR will turn red to indicate charging.
3. The CHARGING LED INDICATOR will turn green when the massager is fully charged.
4. Approximate charging time is 150 minutes.

AREAS OF USE

OPERATING INSTRUCTIONS

1. Check the massager is secure and safe to use.
2. Insert the power plug into a proper wall outlet.
3. Press ON/OFF button to start operate the massage pillow.
4. Press ON/OFF button again to start the automatic heating function. Press ON/OFF button again to switch off the heating function.
5. Hold the ON/OFF button about 3 seconds to switch off the massager.
6. Switch off power and remove power plug from wall outlet.
7. The appliance may be operated for 15 minutes, before a pause of 20 minutes is recommended.

SPECIAL SAFETY PRECAUTIONS

UNSUITABLE USERS

The following users must not use the appliance or should consult their doctor before use:
Anyone undergoing medical treatment or feeling unfit, suffering from hypertension, heart diseases, bran in diseases or thrombosis, malignant turners, pregnant, suffering from osteoporosis, with implanted pacemaker or other medical device.

SAFETY PRECAUTIONS

- Do not operate without carefully reading this operating manual.
- Do not allow children or pets to play around the massager when in use to avoid risk of injury.
- Always wait 20 minutes before reusing the appliance.
- Never use when feeling unwell, if parts are damaged, if smells or smoke occurs during use, during a thunderstorm.
- Do not use while driving.
- Do not use if feeling drowsy and do not fall asleep during use.
- Do not connect to any other power adapter or other electrical equipment at the same time.
- Keep the massager away from water or liquids.
- Do not place the massage heads directly on your hand, elbow, knee or chest.

IMPORTANT: HANDLING OF RECHARGEABLE BATTERIES

- If you come in contact with fluid from the battery cell, rinse the affected areas with water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, open or crush the rechargeable batteries.
- Only use the attached charger.
- Batteries must be fully charged before first time use.
- Correct charging must be observed at all times.
- Fully charge the battery prior to initial use.

DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE)

This symbol on the product or in the instructions means that the electrical and electronic equipment according to the European directive (WEEE) must be disposed of separately from other household waste, as they can be harmful to the environment. The countries in the EU have waste collection systems. Contact your local dealer of authority for more information.

SPECIFICATION

Model no.	13399
Max. operating time	15 minutes
Input voltage	100-240V AC
Output voltage	DC 9V 2A
Frequency	50-60Hz
Power	18W

CONTACT INFORMATION:

F&H Group A/S
GL Skivevej 70
DK-8800 Viborg
DENMARK