# re.

## PROFESSIONAL MASSAGE GUN

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F&H Group A/S  
Gl. Skivevej 70  
DK- 8800 Viborg  
Denmark

## SAFETY INSTRUCTIONS

Please read this instruction manual carefully before first use. It is advisable to keep this instruction manual for future reference.

Use the machine according to the instructions in this manual.

When using the product, all instructions and basic precautions should always be followed, including the following:

**DANGER** - To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning

**WARNING** - To reduce the risk of burns, fire, electric shock, or injury to persons.

1. Use only as instructed. Not for children.
2. Keep massage areas clean and well lit. Cluttered or dark areas invite accidents.
3. Do not operate the massage device in explosive atmospheres, such as in the presence of flammable liquids, gases or dust.
4. Massage devices create sparks which may ignite the dust or fumes.
5. Keep children and bystanders away while operating a massage device. Distractions can cause you to lose control.
6. Stay alert, watch what you are doing and use common sense when operating a massage device. Do not use a massage gun while you are tired or under the influence of drugs, alcohol or medication. A moment of inattention while operating the massage gun may result in serious personal injury.
7. Prevent unintentional starting. Ensure the switch is in the off-position before connecting to the massage battery pack or carrying the device. Carrying massage gun with your finger on the switch or energizing massage gun that has the switch on invites accidents.
8. Do not immerse in water or allow water to enter the ventilation holes of the massage gun.
9. Do not overreach. Keep proper footing and balance at all times. This enables better control of the massage gun in unexpected situations.
10. Only recharge with supplied Massage Gun charger.
11. Dress properly. Do not wear loose clothing or jewellery. Keep your hair, clothing, and gloves away from moving parts. Loose clothes, jewellery or long hair can be caught in moving parts.
12. Carefully examine device and battery before each use.
13. Do not tamper or alter Massage Gun in any way.
14. Do not force the massage gun. Use the correct massage gun for your application. The correct massage device will do the job better and safer at the rate for which it was designed.
15. Do not use the massage gun if the switch does not turn it on and off. Cannot be controlled with the switch is dangerous and must be repaired.
16. Disconnect the battery pack from the massage gun before making any adjustment, changing accessories, or storing the massage gun. Such preventive safety measures reduce the risk of starting the massage gun accidentally.
17. Store the massage gun out of the reach of children and do not allow persons unfamiliar with the massage gun or these instructions to operate the massage gun. A Massage gun is dangerous in the hands of untrained users.
18. Use the massage gun and accessories in accordance with these instructions, taking into account the working conditions and the work to be performed, Use of massage gun for operations different from those intended could result in a hazardous situation.
19. Do not use outdoors.
20. Do not operate where aerosol (spray) products are being used or where oxygen is administered.

21. Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
  22. Do not carry this product by its supply cord or use cord as a handle.
  23. Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. Please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
  24. Use only on a dry and clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region
  25. To avoid pain or discomfort, use the high speed only on the body's soft tissue instead of head or any hard or bony area of the body.
  26. Bruising can occur regardless of speed applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
  27. Keep fingers, hair or any other body part away from the shaft and backside of the head as pinching may occur.
  28. Do not place any objects into the ventilation holes of the massage gun.
  29. Only recharge with supplied 16. BV DC massage gun charger.
  30. Never leave the massage gun operating or charging unattended
  31. This appliance contains batteries that are non-replaceable.
- WARNING:** For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

Adapter: Dong Guan City GangQi Electronic Co., Ltd GQ18-252060-AG

### WARNINGS & CAUTION

1. Do not operate continuously for more than one hour. Allow the device to rest for 30 minutes before re-using.
2. Do not remove screws or attempt to disassemble.
3. To avoid pain or discomfort, use the high speed only on the body's soft tissue instead of head or any hard or bony area of the body.
4. Bruising can occur regardless of speed applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.

**For first time users, please try the product for 10 to 15 minutes first, and then extend massage time if no discomfort occurs. We do not recommend using the product for over 30 minutes per session. Allow device to rest for 30 minutes before re-using.**

### DISCLAIMER

PLEASE DO NOT USE THE MASSAGE GUN OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY.

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis or spondylosis. Recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health.

Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device, but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above.

The deep muscle fascial massage gun is a portable rechargeable vibration device with large torque, long endurance and multi-gear adjustment. It can effectively reduce the pain caused by the accumulation of lactic acid after intense exercise and relieve the discomfort of muscles in the sedentary and long-lying population.


## PACKAGE INCLUDED

The massage gun is a portable rechargeable vibration device. It can effectively reduce the pain caused by the accumulation of lactic acid after intense exercise, and relieve the discomfort of muscles of the sedentary and long-lying population.

Package Included:

1 x Massage gun; 4 x Massage heads; 1 x Charger
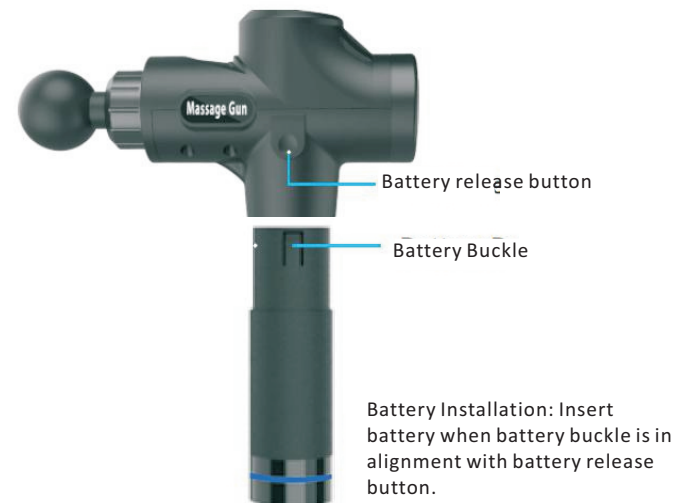
## OPERATING

1. Turn on the power switch at the bottom of the massage gun.
2. Press the ON/OFF button and the Massage Gun starts to run.
3. Press + and - buttons to adjust the speed of the massage gun as your desired.
4. When the battery power is less than 20%, the light line at the bottom of the Massage Gun shows red; when the battery power is less than 10%, the power display window shows E1, it needs to be charged instantly.
5. When the battery is full, the light line at the bottom of the Massage Gun changes from red to green.
6. To remove the battery, lightly press the battery release button with thumb. Firmly pull downwards the battery with the other hand.

## MESSAGE HEAD INTRODUCTION

### BALL HEAD – suits large muscle groups

#### BALL HEAD – suits large muscle groups

The ball massage head is suitable for large muscle groups, like the waist, hips, and legs, as well as muscle relaxation before and after exercise, recommend lightweight fitness enthusiasts who like dancing, yoga, and aerobics.

#### FLAT HEAD – suits most parts of the body

The Flat massage head is suitable for most parts of the body, like chest and large and small muscles. Recommended for fitness people who like to exercise abdominal muscles and pectoral muscles; can be used before and after exercise.

**U-SHAPED HEAD** – suits muscles on both sides of the spine

The U-Shaped Head perfectly avoids the spine and cervical vertebrae of the human body; it is suitable for the shoulders, the cervical vertebrae, and the muscles on both sides of the spine. Not to be used directly on the cervical vertebrae and the spine, which can cause physical damage.

**BULLET HEAD** – suits for meridians

The simple shape on this bullet massage head, which is full of impact experience, it is suitable for the body parts of the meridians, soles, palms, achilles, tendons, back rhomboid muscles; small, but great energy.

#### CHARGING

1. Fully charge the battery up to six hours before first use, either attached or detached from the device.
2. Use supplied adapter to charge the battery
3. Lights flashing means battery is active charging.
4. LED lights stop flashing and remain fully illuminated means the battery is Full charged
5. The battery can be recharged at any time and at any battery level.
6. Have to charge the battery when Led light turns to the red level.
7. The average run-time is + 2 hours depending on speed level and pressure applied during use.

#### MAINTENANCE, CLEANING AND STORAGE

1. Wipe clean with a dry cloth.
2. Store in a supplied carton box in a dry place.
3. Keep the motor ventilation slots clean.
4. Keep all working controls free of dust.
5. Never use water or chemical cleaners to clean your device.
6. For storage or travel, remove battery and store in supplied box.

#### IMPORTANT: handling of rechargeable batteries

- If you come in contact with fluid from the battery cell, rinse the affected areas with water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, open or crush the rechargeable batteries.
- Only use the attached charger.
- Correct charging must be observed at all times.
- Fully charge the battery prior to initial use.

#### DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE)

This symbol on the product or in the instructions means that the electrical and electronic equipment according to the European directive (WEEE) must be disposed of separately from other household waste, as they can be harmful to the environment. The countries in the EU have waste collection systems. Contact your local dealer of authority for more information.

#### SPECIFICATION

Model no.	13401
Power	20W
Motor turn speed	1400-1500-1600-1700-1800-1900-2000-2100-2200-2300-2400-2500-2600-2700-2800-2900-3000-3100-3200-3300/ RPM
Speed	1-20 Speed Setting
Battery Capacity	2600mAh
Rated Voltage	1 x rechargeable 25.2 V
Rated Frequency	20-58Hz
Running Time	2-6 h (depend on speed level)
Charging Time	2.5-4 h
Massage Head	4
Noise Level	50 dB

#### Contact information:

F&H Group A/S  
 Gl. Skivevej 70  
 DK- 8800 Viborg  
 Denmark