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F&H Group A/S
GL Skivevej 70
DK-8800 Viborg
TANSKA

SAFETY INSTRUCTIONS

Please read this instruction manual carefully before first use. It is advisable to keep this instruction manual for future reference.

Use the machine according to the instructions in this manual.

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Use only as instructed. Not for children below 12 years of age.
2. Keep massage areas clean and well lit. Cluttered or dark areas invite accidents.
3. Do not operate the massage device in explosive atmospheres, such as in the presence of flammable liquids, gases or dust.
4. Massage devices create sparks which may ignite dust or fumes.
5. Keep children and bystanders away while operating a massage device. Distractions can cause you to lose control.
6. Stay alert, watch what you are doing and use common sense when operating a massage device. Do not use a massage gun while you are tired or under the influence of drugs, alcohol or medication. A moment of inattention while operating the massage gun may result in serious personal injury.
7. Prevent unintentional starting.
8. Do not immerse in water or allow water to enter the massage gun.
9. Do not overreach. Keep proper footing and balance at all times. This enables better control of the massage gun in unexpected situations.
10. Only recharge with the included supplied charger cable.
11. Dress properly. Do not wear loose clothing or jewellery. Keep your hair, clothing, and gloves away from moving parts. Loose clothes, jewellery or long hair can be caught in moving parts.
12. Carefully examine device and battery before each use.
13. Do not tamper or alter Massage Gun in any way.
14. Do not force the massage gun. Use the correct massage gun for your application. The correct massage device will do the job better and safer at the rate for which it was designed.
15. Store the massage gun out of the reach of children and do not allow persons unfamiliar with the massage gun or these instructions to operate the massage gun. A Massage gun is dangerous in the hands of untrained users.
16. Use the massage gun and accessories in accordance with these instructions, taking into account the working conditions and the work to be performed.

17. Use of massage gun for operations different from those intended could result in a hazardous situation.
18. Do not use outdoors.
19. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
20. Do not carry this product by its supplied cable.
21. Keep the cable away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. Please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
22. Use only on a dry and clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region
23. To avoid pain or discomfort, use the high speed only on the body's soft tissue instead of head or any hard or bony area of the body.
24. Bruising can occur regardless of speed applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
25. Keep fingers, hair or any other body part away from the shaft and backside of the head as pinching may occur.
26. Never leave the massage gun operating or charging unattended.

WARNINGS & CAUTION

1. Do not operate continuously for more than 20 minutes. Allow the device to rest for 30 minutes before re-using.
2. Do not remove screws or attempt to disassemble.
3. To avoid pain or discomfort, use the high speed only on the body's soft tissue instead of head or any hard or bony area of the body.
4. Bruising can occur regardless of speed applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.

For first time users, please try the product for 5 to 10 minutes first, and then extend message time if no discomfort occurs. We do not recommend using the product for over 20 minutes per session. Allow device to rest for 30 minutes before re-using.

DISCLAIMER

PLEASE DO NOT USE THE MASSAGE GUN OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY.

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis or spondylosis. Recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health.


Frail individuals and children should be accompanied by an adult when using any vibration device. These contra indications do not mean that you are not able to use a vibration or exercise device, but we advise you to consult a doctor first.

Ongoing research is being done to the effects of vibration exercise for specific medical disorders. This will likely result in the shortening of the list of contra indications as shown above.

The deep muscle massage gun is a portable rechargeable vibration device with large torque, long endurance and multi-gear adjustment. It can effectively reduce the pain caused by the accumulation of lactic acid after intense exercise and relieve the discomfort of muscles in the sedentary and long-lying population.

PACKAGE INCLUDES:

The massage gun package includes:
1 x Massage Gun, 3 x Massage Heads, 1 x Charging cable and Instruction Manual.

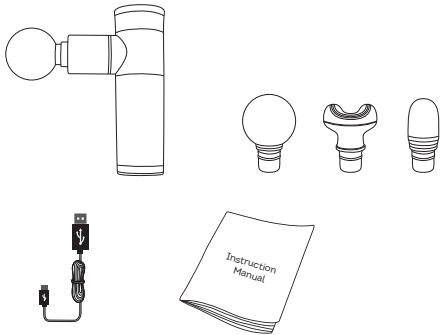OPERATING

Preparation before use


1. Check whether the product and accessories are complete. Do not use them if they are damaged.
2. Check that the battery is sufficiently charged. If the battery is low (red blink), please charge before using.

Charging method

1. Insert the plug on USB cable into the socket of the base, and connect to the adapter USB interface.
 2. The host indicator light is always on pink, indicating charging is in progress. The green light is always on, indicating charging is completed.
- It is normal for the device and power adapter to heat up slightly during charging.

Charging time


1. Charging time is 4-5 hours.
 2. When fully charged it can be used an estimated 8 intervals of 10 minutes before the battery is drained empty. (Due to different massage intensity there will be differences).
 3. When the indicator light flashes red, please charge immediately.
- Unplug the cable from the device and socket when charging is finished.

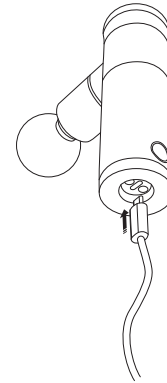HOW TO OPERATE

Schematic diagram of product function configuration

1. Replaceable massage head
2. Light Indicator
3. On/off & Mode switch key
4. Charging interface

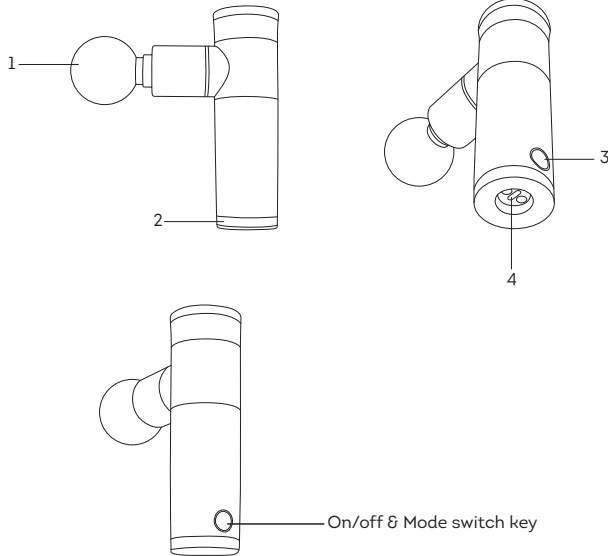HOW TO OPERATE

This product cannot be used while charging. Please unplug power adapter and use it after charging.

1. Install or replace the massage head. Press the power key.
2. On/off machine: Long press the on/off button for two seconds, the machine will start up, the indicator light will turn light green, and the machine will run at a slow speed by default.
3. Intensive adjustment continuous short press on the on/off button to adjust five speed from weak to strong and back to weak in cycles. Different colors indicate different intensity: 1st. level is light green, 2nd. level green, 3rd level is purple, 4th level light pink and 5th level is blue.

- NOTE: The massage gun has a standby function. If pressure isn't applied to the massage head, the speed level stays at 1st level, although the LED light changes colour. The speed level will change once pressure is applied to the massage head.
4. Aim to the parts that need to be massaged. It is important not to use too much force. If the force is too high, the machine will stop working automatically for security protection. After reducing the force or restarting, the machine will work normally.
 5. Automatic shutdown: After 10 minutes of use, the product will automatically shut down.
 6. Shutdown by hand: Long press the on/off button of two seconds.

MESSAGE HEAD INTRODUCTION

BALL HEAD - suits large muscle groups

The ball massage head is suitable for large muscle groups, like the waist, hips, and legs, as well as muscle relaxation before and after exercise, recommend lightweight fitness enthusiasts who like dancing, yoga, and aerobics.

U-SHAPED HEAD - suits muscles on both sides of the spine

The U-Shaped Head perfectly avoids the spine and cervical vertebrae of the human body; it is suitable for the shoulders, the cervical vertebrae, and the muscles on both sides of the spine. Not to be used directly on the cervical vertebrae and the spine, which can cause physical damage.

BULLET HEAD - suits for meridians

The simple shape on this bullet massage head, which is full of impact experience, it is suitable for the body parts of the meridians, soles, palms, achilles, tendons, back rhomboid muscles; small, but great energy.

PRODUCT INSPECTION AND FAULTY OPERATION

Please check if there is a fault in the operation. If the fault cannot be removed or the product still cannot be operated satisfactorily, please turn off and contact your local dealer. Note! Must be repaired by a professional technician. DO NOT disassemble the product yourself!

STORAGE AND MAINTENANCE

Storage

Storage more than 3 months, the machine should be recharged after unpacking and inspection. Do not place it close to the floor, walls or roof. Indoor should be well ventilated, avoid intense sunlight illuminate and corrosive gas erosion.

Transport

The product package is equipped with simple anti-shake, suitable for transportation. Avoid rain, snow splashing, inversion and collision.

Power Inspections

Check the power supply. If the power voltage exceeds the specific range, do not start up and use.

Cleaning

Wipe clean with a dry cloth.

Store in a dry place to prevent moisture, rodents, and pollution.

IMPORTANT: HANDLING OF RECHARGEABLE BATTERIES

- If you come in contact with fluid from the battery cell, rinse the affected areas with water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Risk of explosion! Never leave the massage gun in direct sunlight.
- Do not disassemble, open or crush the rechargeable batteries.
- Only use the attached charger cable.
- Correct charging must be observed at all times.
- Fully charge the battery prior to initial use.

DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT (WEEE)

This symbol on the product or in the instructions means that the electrical and electronic equipment according to the European directive (WEEE) must be disposed of separately from other household waste, as they can be harmful to the environment. The countries in the EU have waste collection systems. Contact your local dealer of authority for more information.

SPECIFICATION

Model no.	24676
Power	10W
Motor turn speed	1528-1885-2290-2495-3650/ RPM
Speed	5 speed intensity
Battery Capacity	2200mAh
Rated Charging Voltage	5V ---
Rated Frequency	230V/50Hz
Running Time	>90 min. (depend on speed level)
Charging Time	4-5 hours
Massage Heads	3
Noise Level	60dB

CONTACT INFORMATION:

F&H Group A/S
GL Skivevej 70
DK-8800 Viborg
DENMARK